

frida baby®

**MediFrida
Accu-Dose Paci**


Quick Start Guide

TAKE A PEEK

**PART TIME PACI.
FULL-TIME MEDICINE
MIRACLE WORKER.**

How to Use:

1. Fill syringe with accurate dose of medicine
2. Twist syringe into pacifier meds channel
3. Place pacifier in baby's mouth with paci plug pointing to either cheek
4. Slowly plunge meds through
5. Twist syringe out & close paci plug so baby can soothe on paci

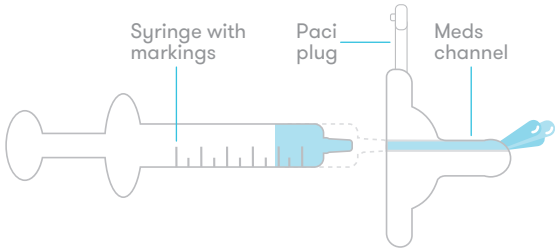

How to Clean:

1. Fill syringe with warm, soapy water
2. Attach pacifier and plunge to clean
3. Rinse + air dry or sterilize

Care + Safety:

Wash pacifier & syringe after each use. Pacifier can be boiled to sanitize between uses. Do not use if syringe markings are worn or unreadable. Follow the medicine bottle's directions for accurate dosage. Do not use if damaged. Do not leave with baby unattended.